

TRIPLE ACTION FOOT MASSAGER



This appliance can be used by children aged from 8 years and above and by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children unless they are supervised.

Please read the user's manual before using the product and keep it for future reference.

A. HEALTH PRECAUTION

Please consult your doctor prior to use if any of the following describes your physical condition

- Pregnant
- Suffer from heart problems, especially if you wear a pacemaker
- Diabetes
- Tuberculosis
- Benign or malignant tumors
- Have concerns about your health

Unless approved by a doctor, do not use massager on areas affected by these conditions if they apply to you :

- hemorrhages,
- wounds,
- swollen areas,
- rash areas,
- burned skin,
- varicose veins,
- open sores,
- poor circulation areas,
- bruises,
- broken skin,
- phlebitis
- thrombosis,
- parts of body affected by paralysis or anesthesia,
unexplained pain in the area where massage is to be applied

within 72 hours following any muscular or joint injury.

B. SAFETY PRECAUTION

READ ALL INSTRUCTIONS BEFORE USING THIS APPLIANCE

1. DANGER

To reduce the risk of electric shock:

- 1) Always unplug this appliance from the electrical outlet immediately when not in use and before putting on or cleaning.
- 2) **DO NOT** use the unit while bathing or sleeping.
- 3) **DO NOT** reach for a product that has fallen into water. Unplug immediately.
- 4) **DO NOT** place in or drop appliance into water or other liquids.
- 5) **DO NOT** place or store appliance where it can fall or be pulled into a tub or sink.
- 6) **DO NOT** stand on this appliance.

2. WARNING


To reduce the risk of fire, electrical shock, burns or injury to persons:

- 1) An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts or cleaning.
- 2) Use heated surface carefully because it may cause serious burns.
- 3) The appliance has a hot surface. Persons insensitive to heat must be careful when using the appliance.
- 4) **DO NOT** use over insensitive skin area or in the presence of poor circulation.
- 5) **DO NOT** operate the appliance if it has a damaged cord or plug, if it is not functioning properly, if it has been damaged, or dropped into water.
- 6) **DO NOT** use outdoors. This appliance is designed for household use only.
- 7) Keep the cord away from heated surfaces, water or liquids.
- 8) **DO NOT** operate under blanket or pillow. Fire, electric shock or injury can occur due to excessive heating.
- 9) **DO NOT** carry the appliance by supply cord or use cord as a handle.
- 10) Use the appliance only for its intended use described in this manual.
- 11) **DO NOT** use attachments not recommended by the manufacturer.
- 12) **DO NOT** drop or insert any object into any opening.
- 13) **DO NOT** operate where aerosol (spray) products are being used or where oxygen is being administered.
- 14) To disconnect, only remove plug from outlet after turning all controls off.
- 15) **CAUTION: DO NOT** use continuously for more than 20 minutes;

If any discomfort is felt, turn off the vibrating foot massager. The next time when you use it, we recommend proceeding with a session lasting up to a maximum of 10 minutes.

- 16) Do not operate the device when it is damaged, or with a damaged cord or plug, or after the product malfunctions, or is dropped or appears damaged in any way. Return the complete product to a service point for inspection, repair or replacement. Always ask for the original spare parts. A similar cable has to be used.

C. DESCRIPTION

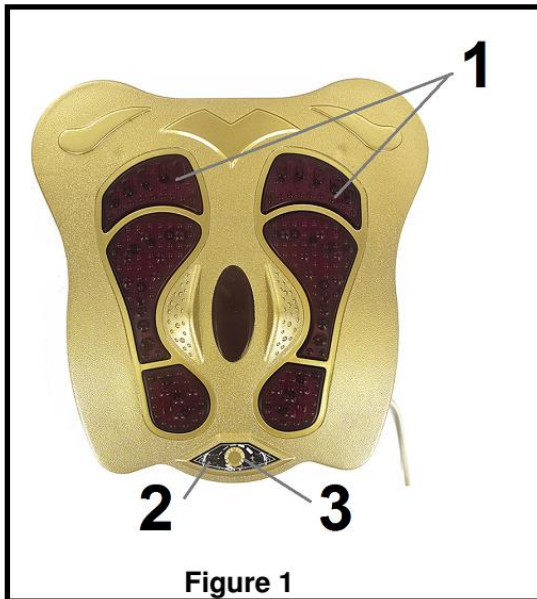

1. Emission area of infrared heat
2. Indicator light
3. ON/ OFF bouton/ Button for setting vibration massage intensity

Figure 1

D. FEATURES

- This innovative foot massager allows 3 actions providing well-being:
 1. Vibrating massage to relax,
 2. Infrared heat to soothe,
 3. Reflexotherapy to stimulate the metabolism.
- Adjustable speed of vibration massage

E. OPERATING INSTRUCTIONS

- 1) Place the unit on the floor firmly, and then plug the cord into the electrical outlet.
- 2) In seated position on a chair or in an armchair or a sofa, place both feet comfortably on the massager.

NOTE: DO NOT stand on massager.

- 3) In order to turn on the massager, turn the button (3) in the clockwise direction: the light indicator and the LEDs are on, and the vibration starts.

Note : The foot massager provides 3 vibration speeds: low, medium, high.

- 4) Then turn the button to adjust the intensity till you get the desired vibration intensity. To reduce intensity, turn the button in the anticlockwise direction.

Note: The foot massager offers several vibration intensities: from low to high. These intensities are obtained by turning the button until you get the desired intensity of vibration.

- 5) After using the foot massager a few minutes, the heat is felt and is continuous until you stop the

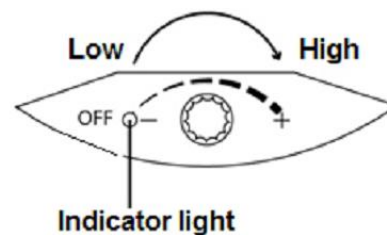


Figure 2

device.

- 6) To stop the foot massager, turn the button anticlockwise until it stops completely. A click is heard to indicate that the foot massager has stopped completely.

F. SPECIFICATIONS

230VAC, 50Hz, 20W

G. MAINTENANCE

- Clean the vibrating foot massager with a dry cloth or towel without water. Its parts that are in direct contact with the skin can be disinfected with alcohol.
- Do not use benzene or any other solvents to clean.
- Store the vibrating foot massager in low temperatures, dry and ventilated place.

H. TROUBLESHOOTING

PROBLEMS	SOLUTIONS
Press ON/OFF, the infrared foot heat emission area does not light up.	<ol style="list-style-type: none">1. First check if the power is ON.2. Turn off power, then start again, and check if it is normal.3. Check if the power plug contact with power socket.

WARNING: You should not dispose of this device with your household waste. A selective collection system for this type of product is implemented by your local authorities. Please contact your local authorities to find out how and where collection takes place. These restrictions apply because electrical and electronic devices contain dangerous substances that have harmful effects on the environment or on human health and must be recycled.



This symbol indicates that electrical and electronic devices are collected selectively. The symbol shows a waste container crossed out with an X symbol. If the power cable is damaged, it must be replaced by the manufacturer, his after sales service, or a person qualified in this area in order to avoid any accidents arising.

If the power cable is damaged, it must be replaced by the manufacturer, his after sales service, or a person qualified in this area in order to avoid any accidents arising.