
i Imported by - Importé par:
GD import (Genius Ideas & Wellys)
190 chemin de la Frayère - 06530 Peymeinade - France
www.gdimport.com
www.genius-ideas.com
www.wellys-care.com

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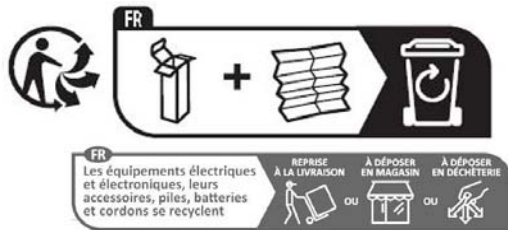GB – TOOTHBRUSH AND SCALER “5 ACCESS”- 028115

Read these operating instructions completely before use and keep them in a safe place for future reference. If the device is passed on to third parties, these operating instructions must also be handed over. If the device is damaged, defective or overheated, it should not be used any further and should be taken to a qualified technician or your dealer. This device is not intended for use by persons (including children under the age of 14) with limited sensory or mental abilities or a lack of experience and / or knowledge, unless they are supervised or received by a person responsible for their safety based on their instructions on how to use the device.

This dental care set is specially designed for the deep cleaning of teeth and for the safe and painless removal of tartar.

The set includes 1 handle, 2 toothbrushes, 2 scaler heads, micro USB cable.

The toothbrush is water resistant (IPX7 standard).

Notes

- At first, there may be slight bleeding from the gums or irritation due to the unusually fast rotating motion. These complaints should only appear for a few days; if they last more than two weeks, stop use immediately and see your dentist.
- If you are undergoing dental treatment, consult your dentist before using the appliance.
- If you have had surgery on your teeth or gums within the past two months, ask your dentist for advice before using this appliance.
- If you wear a pacemaker or other implanted device with built-in electronics or magnets, or if you have other medical concerns, consult your doctor before using the device.

Charging:

To charge, remove the protective cap at the base of the control unit. Connect the toothbrush handle with the USB cable (supplied) to the mains. The green light will flash until charging is complete. Once the full charge is complete the green light will stay on continuously.

Charging time: 1.5 to 2 hours.

Battery life: around 60 days depending on use. When the battery is low the red light will flash after the device is turned off after use.

Note: If the indicator light does not immediately flash green, do not interrupt the charging process. After 10-15 minutes, the indicator light starts flashing.

Warning ! Do not use the device while it is charging!

Use:

The appliance offers three adjustment possibilities :

- Soft ("Gentle") for a gentle and thorough cleaning of sensitive areas
- Strong ("Powerful") for intensive and deep cleaning of all areas
- Comfort ("Whitening") for cleaner, whiter teeth (for occasional use).

1. Choose the desired tip (Brush or Scaler).
 2. Position the tip on the top of the electric handle (rotate the tip so that it snaps into the handle).
 3. Brush tip: Moisten the tip and apply toothpaste. To avoid splashing, put the toothbrush in your mouth before starting the appliance. (It is advisable to change the brush tip every 3 months.)
 4. Scaler tip: After turning on the device, lightly press the cleaning head on the part you want to clean and go back and forth against the scale. If you feel that the tartar is not very hard then it will quickly be powdered. **CAUTION:** Do not use this tip if you have a ceramic prosthesis, it could damage the tip.
- Important: Keep the dental cleaner a sufficient distance from your eyes. Always hold the tip against tooth surfaces when removing tartar. Avoid unnecessary contact of the mouthpiece with soft tissue such as the gums, mucous membranes or skin.
5. Changing modes: When you press the on / off button, your toothbrush automatically switches to "Gentle" mode. To change the mode, successively press the on / off key.
 6. To turn off the appliance, press the on / off button again after the last power mode ("Whitening").

Auto Timer:

The short break at 30 second intervals is meant to remind you when brushing your teeth to evenly clean all four quadrants of the jaw. The device switches off automatically after two minutes.

Note: If you use any other mode for more than 5 seconds, that mode is automatically saved by default and the next time you turn on the device, it will start in that last mode.

Cleaning:

After use, rinse the nozzle thoroughly under water for a few seconds with the handle on. Switch off the handle and remove the mouthpiece. Clean the two parts separately under running water then dry them before reassembling the toothbrush.

Attention:

- Store and charge the device in a safe, dry place that is not in the immediate vicinity of a sink or the like to prevent the device from slipping in water.
- Do not use the product when it is charging.
- The device is intended for private use and not for commercial use. Use the device only as described in the instructions. Any further use is considered inappropriate.
- Do not place the product in water or other liquid for a long time.
- Do not catch a magazine that has fallen in water. Immediately unplug the power plug.
- Do not leave the product on hot surfaces. Do not expose the device to extreme temperatures, large temperature fluctuations, direct sunlight or moisture.
- Do not use additional accessories. Use only the supplied accessories.
- Keep out of the reach of children.
- Do not drop the device or subject it to strong shocks.
- Do not insert any foreign object into any opening of the device. Never attempt to disassemble the device, repair it yourself, or replace any part.
- Do not use the appliance when you are sleepy.
- Do not use this device outdoors or in an environment where sprays or spray materials are used, or where oxygen is administered.
- This product is an individual oral care device and is not intended for multiple patient use in a dental office or clinic.
- Do not put the control unit in contact with alcohol, chlorine, wax, oil or products containing solvents such as perfume, day / sun cream, oil body wash or nail polish remover.
- Do not wrap the cord around the appliance.

Danger Associated with Batteries

- Batteries can be life-threatening when swallowed. Keep batteries and the device out of reach of children. If a battery is swallowed, you must immediately consult medical help.
- Before inserting the batteries, check whether the contacts in the device and on the battery are clean and clean them if necessary. Only use the battery type specified in the technical data.
- Always replace all batteries. Do not use new and used batteries at the same time. Do not use any different battery types, brands or batteries with different capacities. When changing the batteries, check the polarity (+/-).
- Take the batteries out of the unit if they are empty or if you do not use the unit for a prolonged period of time. Thus you will prevent damages that could occur if the batteries leak. If a battery has leaked, ensure that the battery acid does not come into contact with your skin, eyes and mucous membranes. In the event of contact with battery acid, immediately rinse the affected areas with sufficient clean water and immediately consult a physician. Immediately remove a leaked battery from the device. Clean the contacts before you insert new batteries.
- The batteries may not be charged or reactivated with other means, not disassembled, opened, thrown in a fire, submerged in a liquid or short-circuited. Non-rechargeable batteries must not be recharged. The terminals must not be short-circuited.

According to the Electrical and Electronic Equipment Act, you are legally obligated to dispose of old electrical and electronic equipment in an environmentally friendly manner at a municipal disposal facility (recycling center, recycling center, local collection point or return to the dealer). Distributors (dealers) are obliged to take back old devices free of charge. For more information, please contact your local municipal administration. Before disposing of your old electrical or electronic equipment, you are required to remove any batteries or accumulators not enclosed in the equipment, as well as any lamps that can be removed, and to dispose of them safely, separately from the item, at appropriate collection points. The disposal is free of charge.

The following batteries or accumulators are contained in this electrical device:

Battery type	Chemical system
Li-14500-500mAh, 3,7V 1,85W	Lithium

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GB - This product is labelled in accordance with European Directive 2012/19/EU concerning used electrical and electronic appliances (WEEE). Do not dispose of electrical appliances as unsorted municipal waste, use separate collection facilities. Contact your local government for information regarding the collection systems available. If electrical appliances are disposed of in landfills or dumps, hazardous substances can leak into the groundwater and get into the food chain, damaging your health and well-being.

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