
α Alpha Techno GmbH

ELECTRIC MASSAGE CHAIR

AT-8500

USER MANUAL

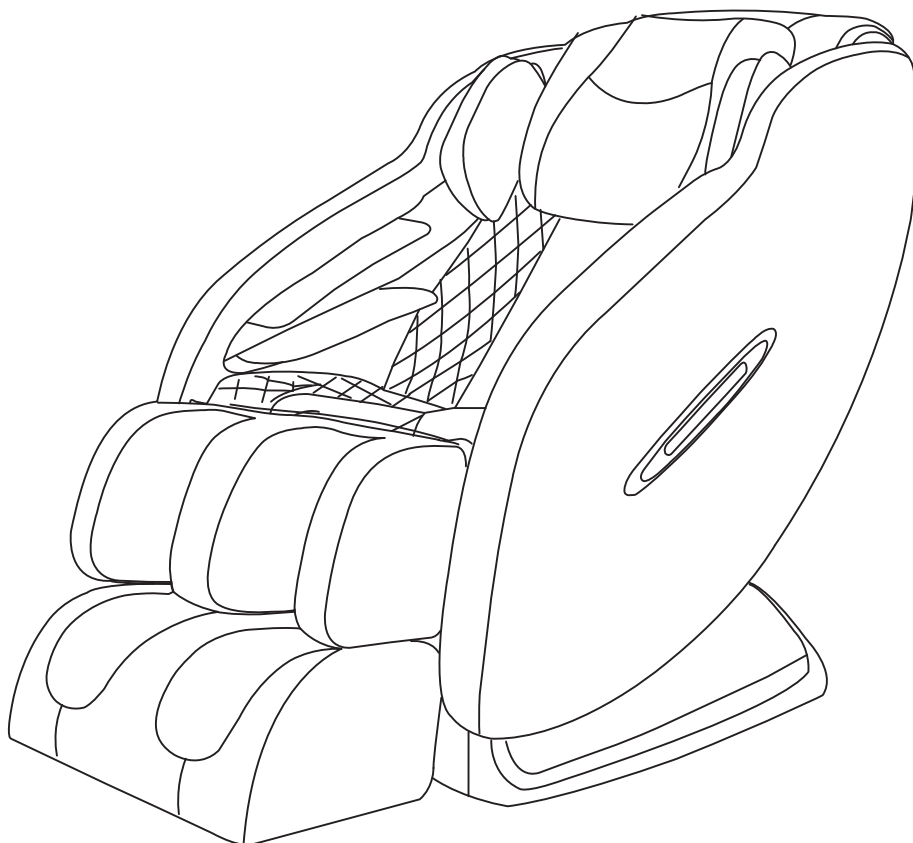Please read this user manual carefully before using the product.
Retain this manual for future reference.

CONTENTS



Thanks for purchasing our massage chair, hopefully the product will supply you really full and relax massage enjoyment. For service issues, please contact your local seller or manufacturer, never dismantle or repair by unauthorized people to avoid any hazard or product damaged.

I	IMPORTANT SAFETY INSTRUCTIONS	1-3
I	PARTS NAME.....	4
I	INSTRUCTIONS ON CONTROL PANEL	5-7
I	HOW TO OPERATE	8
I	CARE AND MAINTENANCE	9
I	SPECIFICATIONS.....	10

IMPORTANT SAFETY INSTRUCTIONS

OPERATING SURROUNDINGS

- Do not use or store the massage chair under damp environment: i.e. bathroom, sauna, spa to avoid any moldy, electric shock, fire or mechanical failure.
- Don't use the massage chair outdoors
- Don't use the massage chair under dusty, smearable and anoxic conditions
- Don't place the massage chair near heater, stove or sunlight radiation to avoid inner components overheating to cause fire
- Make sure to place the massage chair on a flat level ground
- The bracket part should be placed on the plain ground
- To avoid an electric shock, please make sure to connect an earth wire
- Don't use the machine while smoking

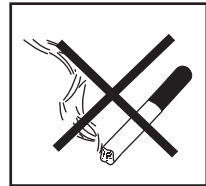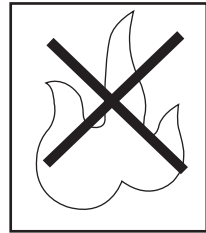FOLLOWING CONDITIONS SHOULD CONSULT DOCTORS BEFORE USE

- Suffered from malignant tumor
- Serious heart disease
- Serious skin disease
- Serious bone hyperplasia
- In pregnancy or during the menses
- Be on the rest cure under the doctor's supervision
- The back bent exceptionally by the accident or disease
- Needing peaceful surroundings
- Arthritis disease
- Implanted pacemaker for medical purpose
- Suffered from bowels disease (i.e. gastritis, enteritis, diarrhea, hepatitis)
- Serious blood circulation trouble
- Consciousness trouble
- Be in hyper pyretic

WARNING

Children under 18 ages or special persons(as psychopath, the mental defects, the paralysis) close are not allowed to operate the massage chair, in case of electric shock, injuries or mechanism malfunctions.

IMPORTANT SAFETY INSTRUCTIONS

For safety reasons and in order to clearly indicate potential dangers should the product be used improperly, instructions have been marked in the following manners and should be strictly noted

⊘ CAUTIONS

Non-professional technicians are not allowed to dismantle or repair to avoid fire, electrical shock or serious hurt

Do not flush to the switch when cleaning to avoid insulation fault, electrical shock or short-circuit

When storing or maintenance, be sure to disconnect the power plug from socket

Do not detach the plug by wet hands

Make sure to shut off and detach the power plug when cleaning or storing product

During cleaning or moving the machine, make sure to detach power cord to avoid fire or electrical shock

Do not place the machine onto an unstable surface to avoid roast, insulation failure, electrical shock or turnover

When storing the machine, make sure to disconnect power plug to avoid fire caused by electricity leaking

Clean power plug periodically to avoid fire due to dust

Do not use it at any unsafe places or on other purposes

It is the main reason of product damage or accident

Do not store the machine onto an unstable surface

Unexpected earthquake or external force will cause damage or injury

Keep supply cord unobstructed at all time, i.e. not twisted, knotted, over flexed or pulled and do not place anything heavy on it

IMPORTANT SAFETY INSTRUCTIONS

CAUTIONS

Do not modify the structure or change components by unauthorized person

Do not dismantle the machine by unauthorized person to avoid electrical shock, short-circuit or fire

Do not store it near a fire or flammable substances

It is one of the reasons of causing product torture, deviation or damage

Should any abnormality occur(i.e. scorching smell), stop using immediately and detach power plug from socket

Keeping operating under an abnormal situation will cause fire or electrical shock, please contact your local seller or authorized service agent for maintenance

Do not place the machine onto an unstable surface to avoid roast, insulation failure, electrical shock or turnover

Only rated voltage must be used

Stop using immediately if supply cord or power plug be damaged or loosen

Be sure to keep children away during operating to avoid any hazard


Do not explore it under sunshine at all time

If the supply cord is damaged,it must be replaced by the manufacturer or its authorized persons for repaire or maintenance to avoid any hazard

PARTS NAME

NAME OF EACH PART

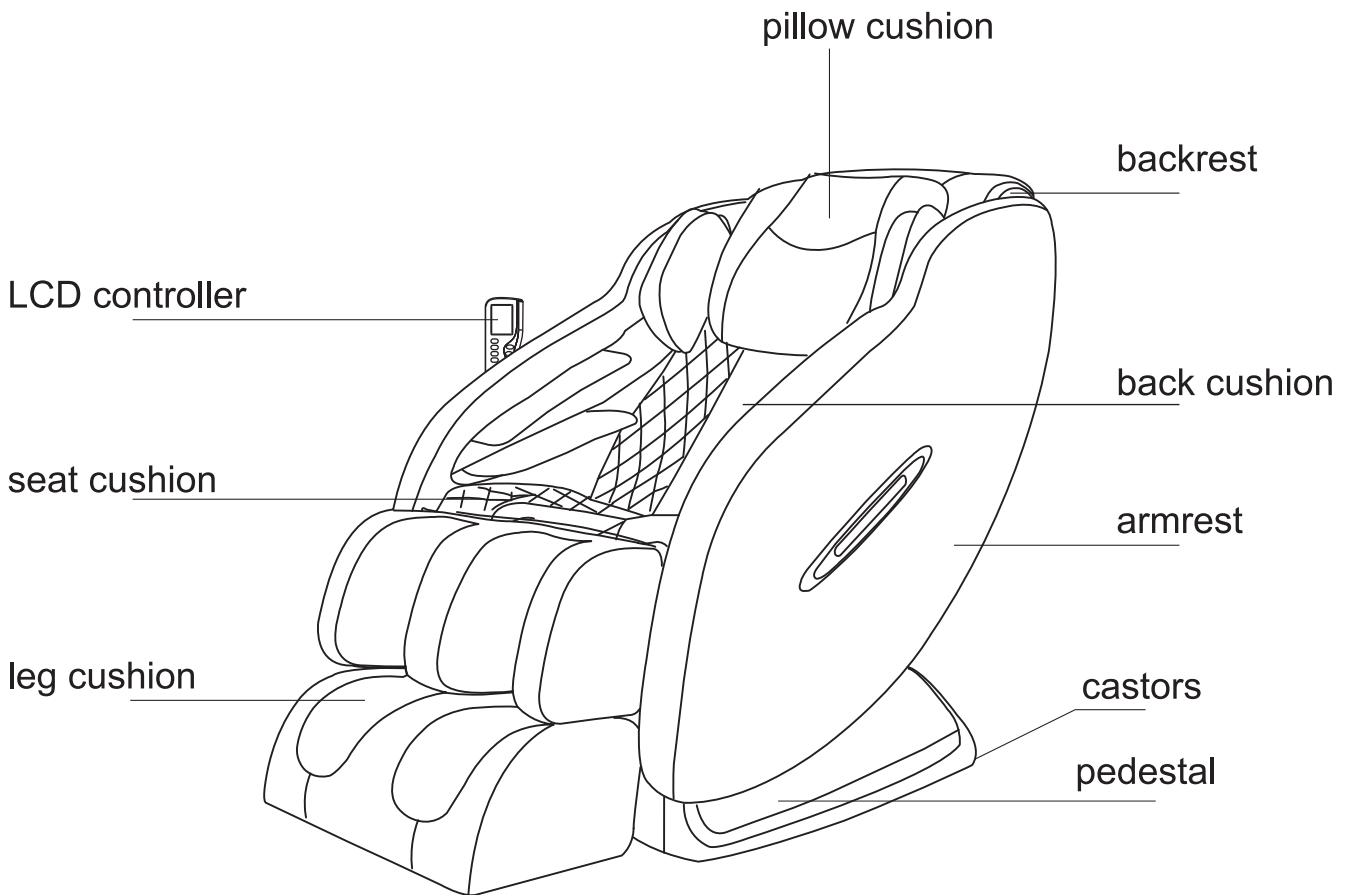WARNING


1. Make sure not to be twisted, knotted or stressed to the power cord or anything heavy on it.
2. The earth wire must be connected properly.
3. Before attaching the plug, make sure all switches must be turned off.
4. If the massage mechanism stays in the center of the backrest, it will reset to the origin position at first.

INSTRUCTIONS ON CONTROL PANEL

HOW TO POWER ON AND POWER OFF

1. Connect the power line and open switch



When connect the power line, press the key to [I] on the box at back of the chair, then the ON/OFF and restore are twinkle, it means the chair has power on.

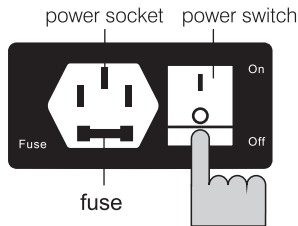2 Sitting the massage chair,press the On/Off key,the chair will start to work

When you sit the massage chair, press the ON/OFF key, the chair will start to work, at this time, the backrest will slow lie and the legrest will slow rise to appoint position, the LCD light,it means the massage chair boot into stand by ON/OFF.

Press the ON/OFF key, it will close all programs,then the backrest and leg rest will return to initial condition.

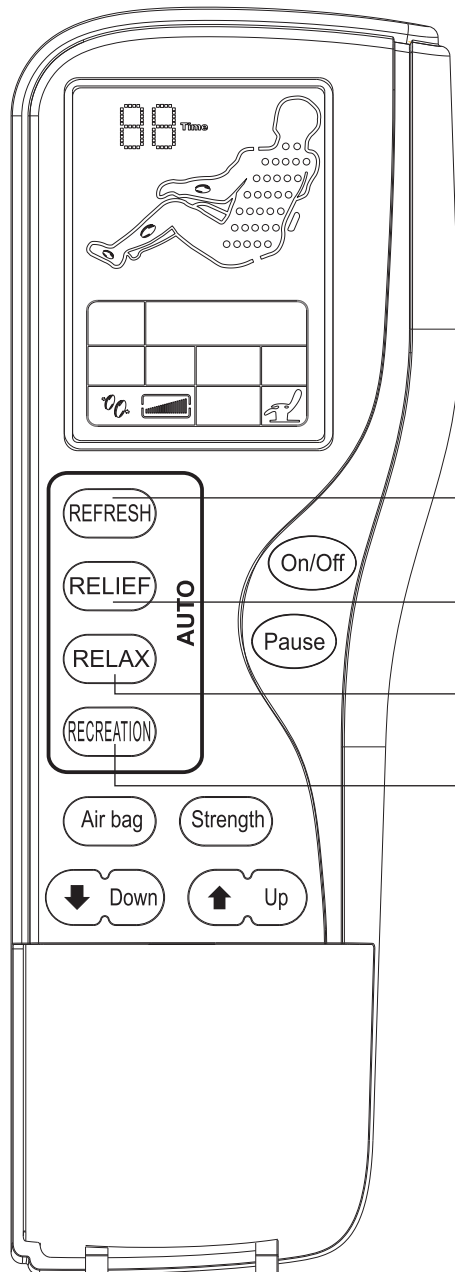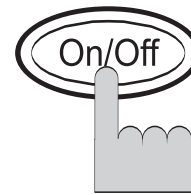Refresh massage program


Relief massage program

Relax massage program

Recreation massage program

Our massage chair has power save setting, LCD will light one minute to automatically extinguish,when you press any key(except On/Off), the LCD will light on.


INSTRUCTIONS ON CONTROL PANEL

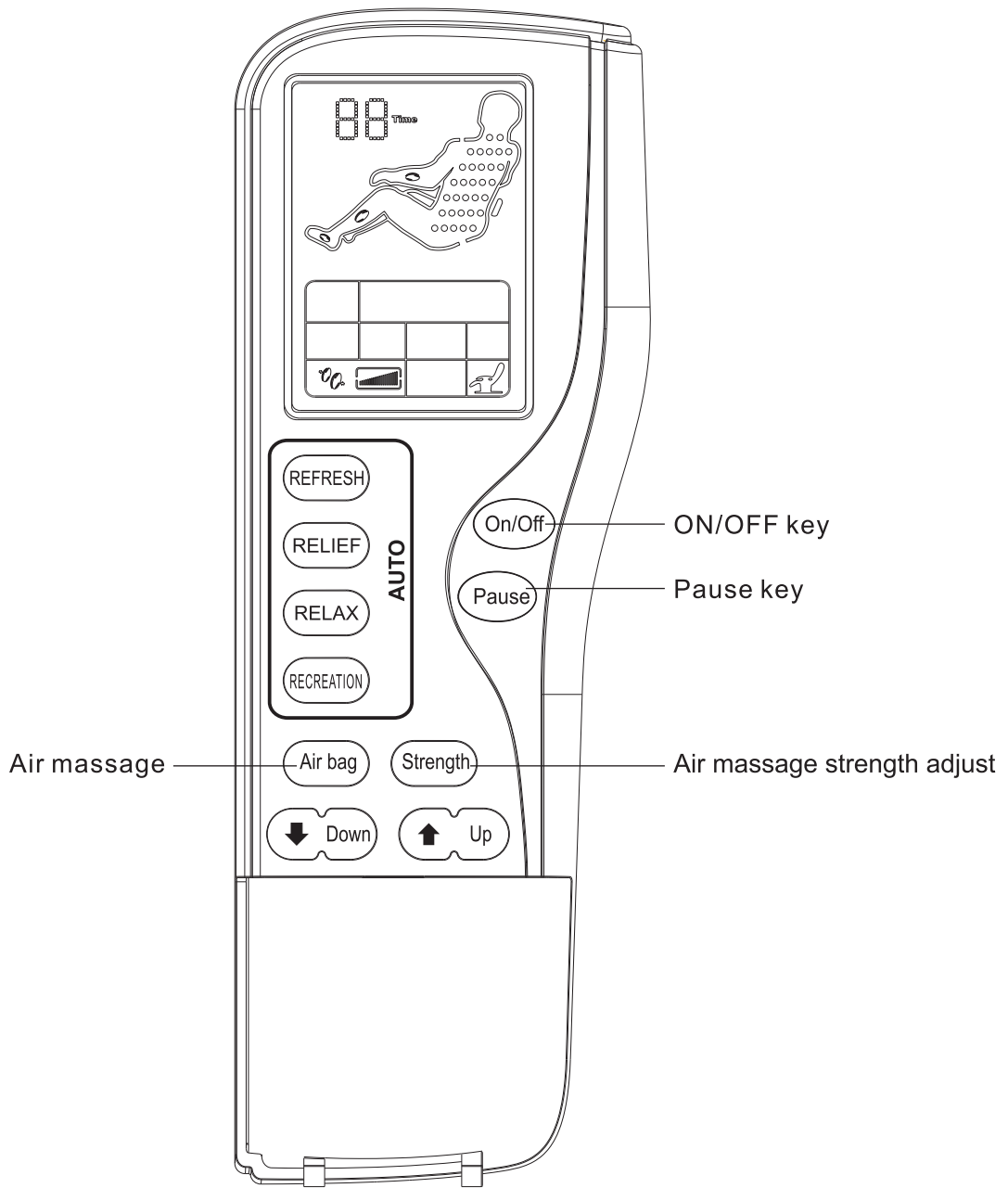3. Adjust the backrest and legrest


Press backrest and legrest key readjust the angle you feel comfortable. When you press backrest up/down and legrest up/down ,they are will laid or risen ,release the press will stop to adjust the angle.

Backrest & Legrest

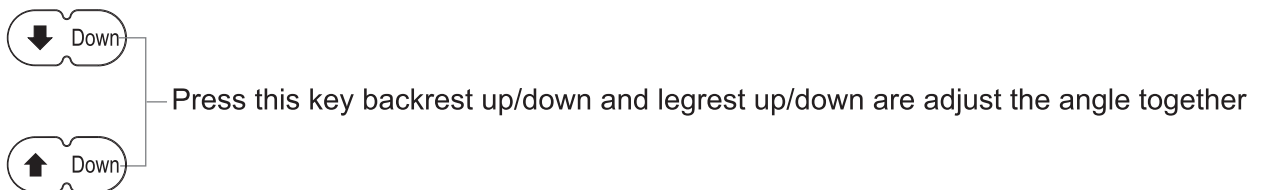INSTRUCTIONS ON CONTROL PANEL

INSTRUCTIONS OF BUTTONS

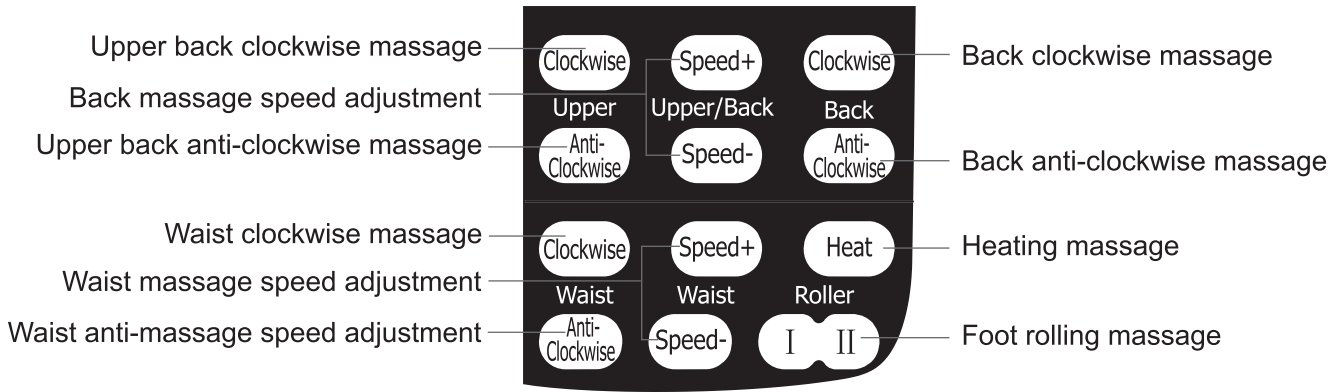4. Stop or finish massage

Normal stop : the timing function automatically fix time is 15 minutes,time out,all the program will stop and keep the setting posture.

Force stop: when press the pause during massage, it will stop massage immediately, press again to continue.

Avoid the legrest damage,please lie the legrest to the initial condition when you stop massage.

5. Start massage

Auto massage program

REFRESH, RELIEF ,RELAX AND RECREATION

Auto massage program has different massage technique:

Kneading and Air press massage.

KEYS	DESCRIPTION	MAIN FUNCTION
REFRESH	Refresh massage program	kneading massage, air massage and tapping massage. It can adjust the angle during massage.
RELIEF	Relief massage program	Back massage and leg massage. It can adjust the angle during massage.
RELAX	Relax massage program	Slow massage and air massage. It can adjust the angle during massage.
RECREATION	Recreation massage program	Back massage and leg massage. It can adjust the angle during massage.


All of the auto massage program, the timing is 15 minutes ,air bag strength is medium,also you can press the STRENGTH to adjust the strength of air bag.

All the auto massage program is stand alone program.

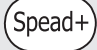HOW TO OPERATE

INSTRUCTIONS OF PARTS AUTO MASSAGE PROGRAM

CLOCKWISE AND ANTI-CLOCKWISE

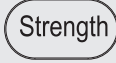
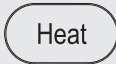
KEYS	DESCRIPTION	MAIN FUNCTION
 Upper 	Upper back clockwise massage	Massage speed adjustable
	Upper back anti-clockwise massage	Massage speed adjustable
 Upper+Back 	Upper back or whole back massage speed increase	3 levels massage speed can be choosed
	Upper back or whole back massage speed decrease	3 levels massage speed can be choosed
 Back 	Whole back clockwise massage	Massage speed adjustable
	Whole back anti-clockwise massage	Massage speed adjustable
 Waist 	Waist forward massage	Massage speed adjustable
	Waist reverse massage	Massage speed adjustable
 Waist 	Waist massage speed increase	3 levels massage speed can be choosed
	Waist massage speed decrease	3 levels massage speed can be choosed

The timing of auto massage program is 15 minutes,the speed is medium,and you can press SPEED to adjust speed.

HOW TO OPERATE

FUNCTION MESSAGE

AIR BAG AND TIMER

KEYS	DESCRIPTION	MAIN FUNCTION
	Air bag massage	It' s used for on/off air bag massage.
	The adjust key of air bag strength	Choose five levels of air bag strength. It should fit air bag massage program.
	Foot rolling massage key	Choose two levels of foot rolling massage.
	Waist heating massage	It' s used for on/off waist heating massage.

All the above is stand alone massage,it can open and close in any massage program.

CARE AND MAINTENANCE

STORAGE

- Do not remove the backrest form the chair ,touch the internal spare parts of massage chair may reduce electric shock.
- Do not remove the internal spare parts or repair the massage chair by yourself .It may reduce electric shock and risk of fire.
- Please pull out the power line when you are maintenance , if not it may reduce electric shock and hurt.
- The massage chair should keep clean all the time , and set in the dry place.
- If the massage chair not use for a long time ,please cover cloth to avoid dust.
- Avoid direct sunlight.

THE MAINTENANCE OF MESSAGE CHAIR

- When you clean the massage chair please use neutral detergent and special detergent for clothing
- If the backrest cushion , pillow cushion and legrest cushion are dirty ,please use neutral detergent and air seasoning
- The plastic,piping and leather could use detergent ,then clean with dry cloth
- Do not use abrasive cleaning agents like benzene or thinner to clean the controller, this may result in malfunction
- Do not spray the insecticide to the massage chair,it may reduce electric shock and risk of fire
- Please use dry cloth cleaning to avoid trouble

THE MAINTENANCE OF CONTROLLER

- Avoid source of heat and source of ignition
- Keep dry and avoid electric circuit
- Please cleaning by soft and dry cloth,and dust elimination by small banister brush
- Prohibit the use of various volatile solvent(i.e.tolueue ,oil,lacquer thinner etc.)

CARE AND MAINTENANCE

TROUBLESHOOTING

Please surveied the massage chair when you are used it often

Please checked up the massage chair when you no used for a long time

SITUATIONS	SOLUTION
The massage chair doesn' t work ,ON/OFF and Restore key don' t light up	1.Check whether the plug is well connected and working well 2.Check whether the switch on the massage chair is turnen to [I]positon
Any sound during work	Under those sound and feeling is normal phenomenon 1. The sound when tapping and kneading 2. The different after adjust “speed” 3. Internal sound when massage 4. Massage wheel grated on the cloth 5. Expansive sound and exhaust sound during air massage 6. The sound of adjust legrest 7. The sound of motor when exert the wight of load

EMERGENCY TREATMENT

1. An exception occurred during the massage chair work,please close the switch and pull out the power line.
2. If the power line is broken,please notify our customer service to processing.

Warning: Except professional you don' t dismantling,repair and change the massage chair,it may reduce electric shock and risk of fire.


SPECIFICATIONS

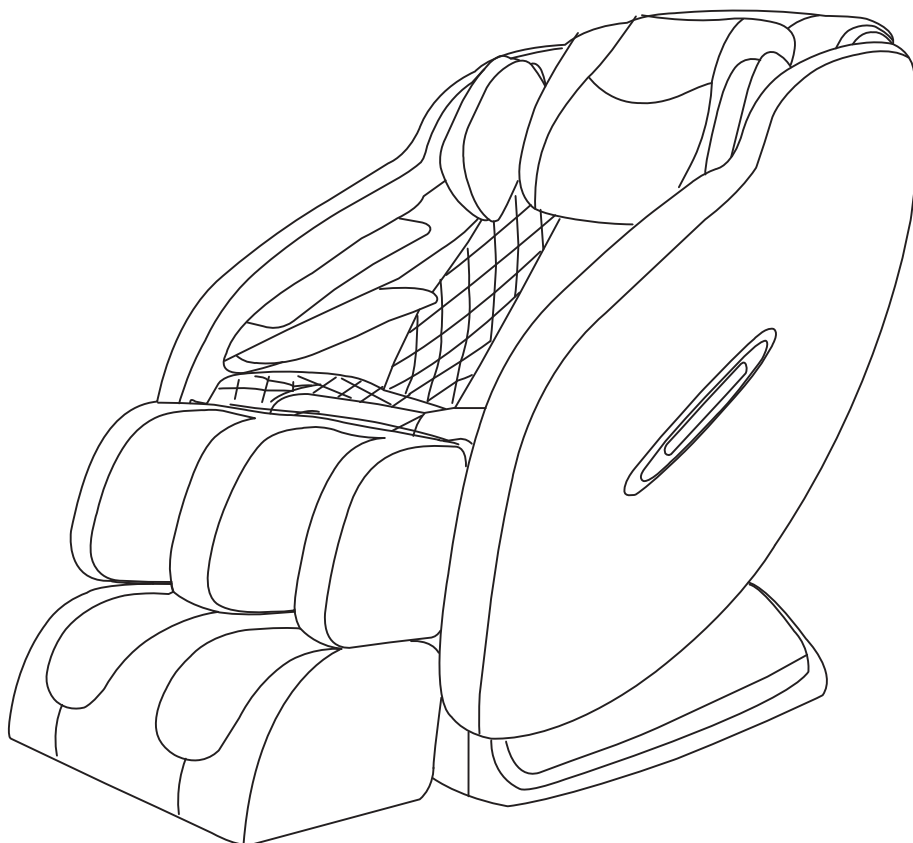Description	Massage Chair
Model	AT-8500
Rated supply	230V~
Rated Frequency	50Hz
Rated Power	100W
Rated working hour	15 min

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